

Improved Design

LoFric[®] Hydro-Kit[™] is small and convenient with proven performance. It's the all-in-one catheter system that's ready when you are, wherever you are!

16" length



Integrated sterile water sachet.
Fold and press to release the water and activate the coating.

Now, a more convenient bag.
Longer and narrower 1000mL urine collection bag.

Convenient loop.
Use the loop to hang the bag or insert your thumb to hold it.


Now, easier to empty.
An improved tear point makes it easier to open and empty the bag.

No-touch feature.
The handling aid provides a firm grip without touching the catheter.

Low-friction catheter with Urotonic[™] Surface Technology.
Minimizes the risk of infections and maximizes comfort.


Now, less waste and easier to handle.
The outer sleeve has been removed and two welded points keep the kit folded until ready for use.

8" length

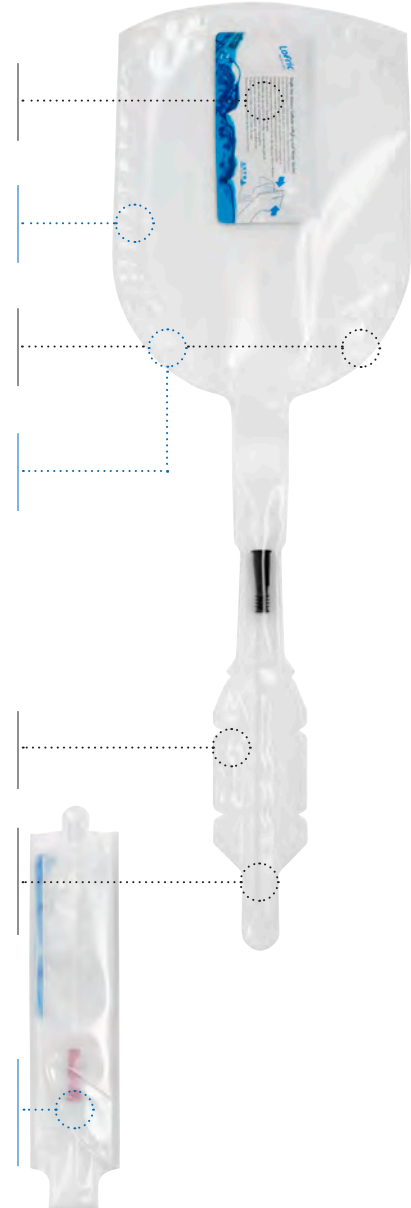


Not made with natural rubber latex.


Instructions

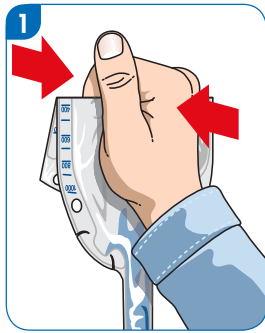

1. Unfold product and hold upright. Fold the water sachet and squeeze.

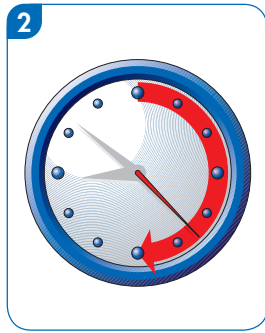

2. Let the water run down to the catheter and soak for 30 seconds. The loop can be used to hang up the product.

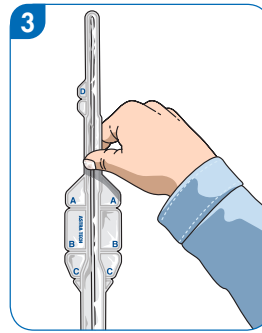


3. Turn the product upside down to allow the water to drain into the collection bag.



4. Tear off at "A" to uncover the tip of the catheter. Before removing it, use this piece to help pull the catheter out until the funnel end has a tight seal against the neck of the bag. This will prevent leakage.



5. Tear off at "B". Use section "B" as a "no-touch" handling aid. Catheterize.

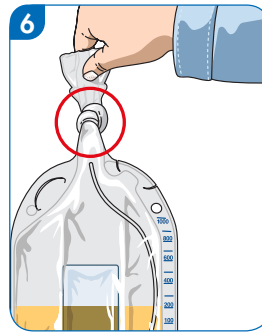

6. To discard after use tear off at "C" and discard the catheter then tie a knot to seal the bag.



7. Tear at the indentation to empty the collection bag.

Product Range

Item #	Type	Length	Size	Color Code
9830840	LoFric Hydro-Kit, Male	16"	8FR	Blue
9831040	LoFric Hydro-Kit, Male	16"	10FR	Black
9831240	LoFric Hydro-Kit, Male	16"	12FR	White
9831440	LoFric Hydro-Kit, Male	16"	14FR	Green
9831640	LoFric Hydro-Kit, Male	16"	16FR	Orange
9831840	LoFric Hydro-Kit, Male	16"	18FR	Red
9871040	LoFric Hydro-Kit, Coudé	16"	10FR	Black
9871240	LoFric Hydro-Kit, Coudé	16"	12FR	White
9871440	LoFric Hydro-Kit, Coudé	16"	14FR	Green
9871640	LoFric Hydro-Kit, Coudé	16"	16FR	Orange
9871840	LoFric Hydro-Kit, Coudé	16"	18FR	Red
4230840	LoFric Hydro-Kit, Female	8"	8FR	Blue
4231040	LoFric Hydro-Kit, Female	8"	10FR	Black
4231240	LoFric Hydro-Kit, Female	8"	12FR	White
4231440	LoFric Hydro-Kit, Female	8"	14FR	Green
4231640	LoFric Hydro-Kit, Female	8"	16FR	Orange
4231840	LoFric Hydro-Kit, Female	8"	18FR	Red
4210840	LoFric Hydro-Kit, Pediatric	8"	8FR	Blue
4211040	LoFric Hydro-Kit, Pediatric	8"	10FR	Black

Intended use

Catheter and collection bag for intermittent catheterization and bladder emptying.

Material composition

Hydrophilic surface coating: poly(vinyl pyrrolidone)/PVP and table salt/NaCl

Catheter: 16" products: medical grade poly(vinyl chloride)/PVC
8" products: polyolefin-based elastomer/POBE

Plastic packaging: poly(ethylene butylacrylate)/EBA

Water pocket: Laminate foil of poly(ethylene terephthalate)/PET, polyethylene/PE and Aluminium, Sterile water

Sterilization and shelf life

All catheters are sterilized by irradiation. Information about the recommended best-before date is printed on the package.

Package

20 catheters per carton and 2 cartons per shipping box.