


LoFric Origo Pro Coudé – quick guide

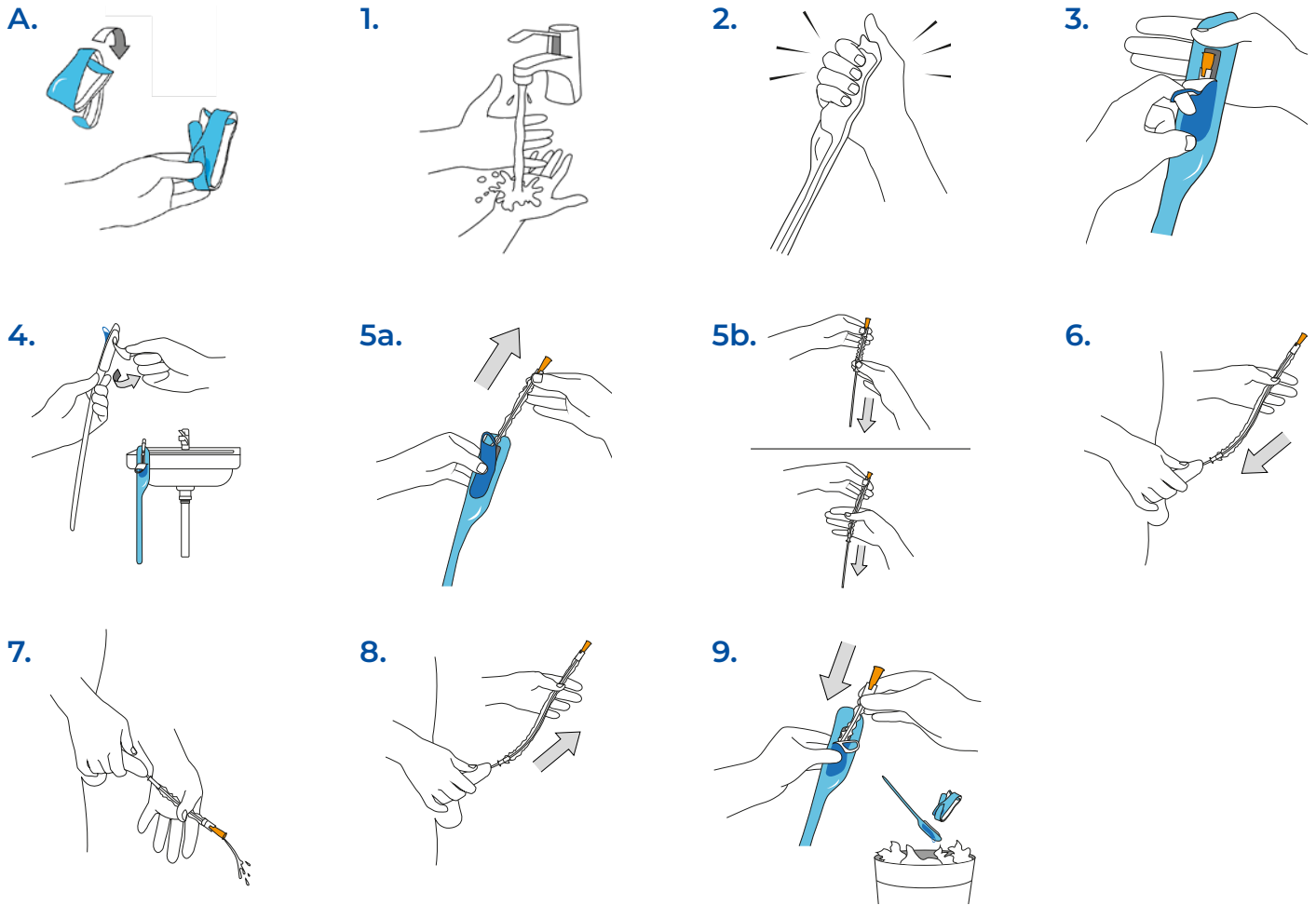

A. Tips for storing before use. Foldable to pocket size.
Use the slit to keep the catheter folded.

1. Wash your hands thoroughly with soap and water.
2. Squeeze to release the salt solution and the catheter is ready to use.
3. Use the loop to pull down and open.
4. Optional: Use the adhesive tab on the reverse side to attach the product to a dry, clean surface.
5. A. After opening, gently pinch and hold the sleeve grip to automatically pull down the protective sleeve when taking out the catheter. Use the protective sleeve to control insertion without having to touch the catheter tube.

B. OPTIONAL: Take out the catheter. Pull the sleeve grip down to cover the catheter tube. Use the protective sleeve to control insertion without having to touch the catheter tube.

6. Lift the penis towards the stomach to straighten the urethra. Slowly insert the catheter into the urethra until urine starts to flow, insert slightly more to ensure a steady stream.
7. Angle the penis down as urine begins to flow through the catheter.
8. Wait until the urine flow stops, to ensure complete bladder emptying. Slowly withdraw the catheter.
9. Put the catheter back in the package and dispose appropriately (local regulations may vary).

IMPORTANT INFORMATION

Use only the salt solution included to wet the catheter. Do not place the catheter on paper or cloth as these materials can absorb the wet coating.

Important information

LoFric® Origo™ Pro Coudé, single use, hydrophilic urinary catheter, with salt solution for activation and multiple eyelets, designed to allow for complete bladder emptying without the need for repositioning.

For short and long term bladder management with intermittent urinary catheterization.

For adults, adolescents.

LoFric catheters are for prescription use only.

Ensure to follow the advice and training on the catheterization procedure given by your Healthcare Professional, before use.

Please refer to instructions for use included in the package, for important information and complete instructions.

Contact your prescriber if you experience difficulties.

At Wellspect, we develop innovative products and solutions that make a real difference in the lives of people managing continence care needs, empowering them to live each day to the fullest. We have been a trusted leader in healthcare for over 40 years with our product brands LoFric®, Navina™ and Surity™, continuously expanding access to our user-friendly solutions. Every day, we strive to maintain the highest quality standards while minimizing the environmental impact of our products. By working closely with healthcare professionals and people who rely on our solutions, we aim to improve clinical outcomes and enhance quality of life for people around the world.

Wellspect. A Real Difference.

For more information about our products and services, please visit [Wellspect.us](https://wellspect.us).

Join the conversation on LinkedIn, X, Facebook and Instagram.

wellspect.us

Wellspect HealthCare, 1235 Friendship Rd, Suite 205, Braselton, GA 30517. Toll-free: +1 877 456 3742 .

LoFric



Manufacturer: Wellspect HealthCare, Aminogatan 1, P.O. Box 14, SE-431 21 Mölndal, Sweden.
Phone: +46 31 376 40 00. www.wellspect.com

CE 2797

R_x Only